



FRONTIER™ 5706P

Multi Centrifuges



Reliable and Compact Multi-Function Low-Speed Centrifuge for General Purpose

Frontier 5706P Multi Centrifuge is a small benchtop multi-function centrifuge ideal for a multitude of separation applications. With versatile standard rotor options, the FC5706P seamlessly processes both round and conical bottom tubes, including common blood and urine tubes. Enhanced by user-friendly design, this centrifuge features an intuitive interface that streamlines operation and maximizes efficiency, all while conserving valuable bench space.

Features:

- Designed for ease of operation, Frontier 5706 Multi Centrifuge features convenient control on an intuitive interface, which allows for responsive operation – even when wearing gloves.
- Enhanced adaptability to satisfy different customers' requirements in general laboratory applications. The FC5706P can hold up a variety of different tubes, starting from 50 ml to smaller volume 1.5 ml tubes. The Frontier 5706P offers 5 rotor options including a fixed angle rotor suitable for 18×1.5/2.0 ml tubes and a swing-out rotor that can hold 4×10 ml blood and urine tubes, boosting its versatility for various laboratory needs.
- Enclosing all its functionalities within a small footprint, the FC5706P saves valuable space on the lab bench.

FRONTIER™ 5706P *Multi Centrifuges*

Designed for a User-Friendly Experience

- Easy-to-use touch-wheel interface for parameter adjustments
- Convenient interface shows both set and actual parameters
- One-touch toggling between RPM/RCF and Acc/Dec
- Double-line LCD display with light blue back light
- Audible signal at the end of each run with a melody
- 10 acceleration/deceleration settings
- 99 memory sets to quickly store and recall commonly used parameters



Safety Features Enclosed Within Small Footprint

- Microprocessors control to help prevent accidents and protect both the user and the equipment
- Rotor imbalance sensor system stops the unit when it detects an imbalanced load
- Sturdy mechanical lid with locking system
- LED status to indicate the operational condition of the centrifuge like run completed (green), standby (yellow) or centrifugation (red)



Flexible Rotor Selections and Kits to Cover a Wide Range of General Laboratory and Diagnostic Applications

- Maintenance-free DC motor produces a maximum G-force of $4427 \times g$
- Speeds up to 6000 rpm and quick spin for short runs
- Additional fixed angle rotor (18 × 1.5/2.0 ml microtubes) and swing out rotor (4 × 10 ml blood tubes) for enhanced flexibility



FRONTIER™ 5706P Multi Centrifuges

Specifications

Model		FC5706P
Speed Range		200 – 6000 rpm; 10 rpm/set
Maximum RCF		4427 × g; 10 g/set
Maximum Capacity (Rotor)		Angle Rotor: 6 × 50 ml
Temperature range		Air cooled
Refrigerant		N/A
Running Time		Quick Spin, 10 sec to 99 hr 59 min, or continuous
Noise level (depending on the rotor)		≤ 56 ±2 dB(A)
Acc/Dec Levels		10/10
Allowable Density at Maximum Speed		1.2 g/ml
Allowable Kinetic Energy		2427 Nm
Rated voltage		100-240 V/50-60 Hz
Voltage Fluctuation		± 10 %
Rated current		2.5-1.0 A
Power continues		70 W
Dimensions (W×D×H)		323 x 427 x 271 mm / 12.7 x 16.8 x 10.7 in
Net Weight (without rotor)		14.2 kg / 31.3 lb
Shipping Dimensions (W×D×H)		370 x 540 x 365 mm / 14.6 x 21.2 x 14.4 in
Shipping Weight (without rotor)		17.6 kg / 38.8 lb
Ambient conditions (EN/IEC 61010-1)	Environment	For indoor use only
	Altitude	Use up to an altitude of 2000 m
	Ambient temperature	5 °C up to 35 °C
	Max. relative humidity	Max. relative humidity 80 % for temperatures up to 31 °C Decreasing linearly to 50 % relative humidity up to 35 °C.
EMC		EN/IEC 61326-1 Class B emissions, Basic immunity FCC Class B emissions
Not suitable for use in hazardous environments.		

FRONTIER™ 5706P Multi Centrifuges

Accessories

	Order No.	30130877
	Description	Rotor Angle 12 × 15 ml RB/FA
	Included in Kit	31012113
	Max. Speed	6000 rpm
	Max. RCF	4427 × g

Accessories

Order No.	Description	Vol./Tube	Diameter	Tube Type	No.of Tubes/ Rotor
Without Adapter	Rotor Angle 12 × 15 ml RB/FA	15 ml	17 mm	Round Bottom/Conical Bottom	12
30211180	Adapter 1 × 6-10 ml D16 mm Vac/Sar, 2 pk	6-10 ml	16 mm	Vacutainer®/Sarstedt	12
30130889	Adapter 1 × 4-7 ml D13 mm Vac/Sar, 2 pk	4-7 ml	13 mm	Round Bottom/Vacutainer®/Sarstedt	12
30130890	Adapter 1 × 1.8-5 ml D13 mm Vac/Sar, 2 pk	1.8-5 ml	13 mm	Round Bottom/Vacutainer®/Sarstedt	12
30130886	Adapter 1 × 1.5/2.0 ml D11mm, 6 pk	1.5/2.0 ml	11 mm	Microtube	12

	Order No.	30130878
	Description	Rotor Angle 6 × 50 ml RB/FA
	Included in Kit	31012115
	Max. Speed	6000 rpm
	Max. RCF	4427 × g

Accessories

Order No.	Description	Vol./Tube	Diameter	Tube Type	No.of Tubes/ Rotor
Without Adapter	Rotor Angle 6 × 50 ml RB/FA	50 ml	29 mm	Round Bottom/Conical Bottom	6
30130891	Adapter 1 × 30 ml D26 mm, 2 pk	30 ml	26 mm	Round Bottom	6
30130892	Adapter 1 × 16 ml D18 mm, 2 pk	16 ml	18 mm	Round Bottom	6
30130893	Adapter 1 × 15 ml D17 mm, 2 pk	15 ml	17 mm	Round Bottom	6
30130894	Adapter 1 × 15 ml D17 mm RB/FA, 2 pk	15 ml	17 mm	Round Bottom/Conical Bottom	6

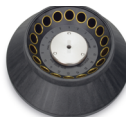
	Order No.	30130880
	Description	Rotor Swing out 6 × 5 ml RB
	Max. Speed	4 000 rpm
	Max. RCF	1 878 × g

Accessories

Order No.	Description	Vol./Tube	Diameter	Tube Type	No.of Tubes/ Rotor
30130895	Buckets For 30130880, 2 pk*	5 ml	13 mm	Round Bottom/Vacutainer®/Sarstedt	6

FRONTIER™ 5706P Multi Centrifuges

	Order No.	31013397				
	Description	Rotor Swing Out 4 × 10 ml				
	Included in Kit	31012114				
	Max. Speed	4200 rpm				
	Max. RCF	2406 × g				
Accessories						
Order No.	Description	Vol./Tube	Diameter	Tube Type		No.of Tubes/ Rotor
31021231	Rack 1 × 1.8-10 ml D16 mm Vac/Sar, 4 pk**	1.8-10 ml	16 mm	Round Bottom/ Vacutainer®/ Sarstedt		4

	Order No.	83041499			
	Description	Rotor Angle 18×1.5/2.0 ml			
	Max. Speed	6000 rpm			
	Max. RCF	2938×g			
Accessories					
Order No.	Description	Vol./Tube	Diameter	Tube type	No.of tubes / rotor
Without Adapter	Rotor Angle 18×1.5/2.0 ml	1.5/2.0 ml	11 mm	Microtube	18
30130884	Adapter 1×0.2/0.4 ml D6 mm, 6 pk	0.2/0.4 ml	6 mm	Microtube	18
30130885	Adapter 1×0.5 ml D8 mm, 6 pk	0.5 ml	8 mm	Microtube	18

*Buckets included

**Racks included

Shortcut list:

Vac – Vacutainer® Sar – Sarstedt FA – Conical bottom / Falcon® FB – Flat bottom RB – Round bottom

FRONTIER™ 5706P *Multi Centrifuges*

Other Standard Features and Equipment

Microprocessor control, Imbalance switch-off, Audible signal at the end of each run with melody options, Emergency release, 99 memory sets, Ten acceleration and ten deceleration rates, Pre-selection of running time from 10 sec to 99 hr 59 min or continuous, Quick spin for short runs.

Compliance

See Declaration of Conformity available CE Declarations | OHAUS

Outline Dimensions

