

Instruction Manual

HI 991001

Extended Range Waterproof pH / Temperature Meter



WARRANTY

All Hanna Instruments meters are warranted for two years against defects in workmanship and materials when used for their intended purpose and maintained according to instructions. The probes are warranted for a period of six months.

This warranty is limited to repair or replacement free of charge. Damage due to accidents, misuse, tampering or lack of prescribed maintenance are not covered.

If service is required, contact the dealer from whom you purchased the instrument. If under warranty, report the model number, date of purchase, serial number and the nature of the problem. If the repair is not covered by the warranty, you will be notified of the charges incurred. If the instrument is to be returned to Hanna Instruments, first obtain a Returned Goods Authorization Number from the Customer Service department and then send it with shipment costs prepaid. When shipping any instrument, make sure it is properly packaged for complete protection.

All rights are reserved. Reproduction in whole or in part is prohibited without the written consent of the copyright owner.

Dear Customer,

Thank you for choosing a Hanna Instruments Product.

Please read this instruction manual carefully before using the instrument. If you need additional technical information, do not hesitate to e-mail us at tech@hannainst.com.

PRELIMINARY EXAMINATION

Remove the instrument from the packing material and examine it carefully to make sure that no damage has occurred during shipment. If noticeable damage is evident, notify your dealer.

Note: Save all packing material until you are sure that the instrument functions correctly. All defective items must be returned in the original packing together with the supplied accessories.

GENERAL DESCRIPTION

HI 991001 is a portable, microprocessor-based pH and temperature meter.


Main features include: extended pH and temperature ranges; waterproof and compact casing; large dual-line display; low battery detection; automatic pH calibration at one or two points within two memorized buffer sets (standard or NIST); selectable temperature unit (°C or °F).

The HI 1296D pH electrode has a built-in temperature sensor for simultaneous temperature compensated pH and temperature readings, and also contains an integral pH sensor preamplifier to provide measurements impervious to noise and electrical interferences.

Each meter is supplied with:

- HI 1296D pH/temperature probe with titanium body and DIN connector and 1 m (3.3') cable
- Calibration & cleaning solution sachets
- 3 x 1.5V AAA batteries
- Instruction manual.

LCD DESCRIPTION

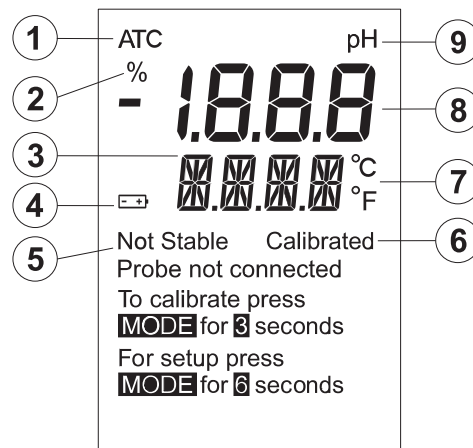

1. Automatic Temperature Compensation indicator
2. Battery percentage (visible at power up)
3. Secondary display
4. Low battery indicator
5. Stability indicator
6. pH calibration indicator
7. Selectable temperature unit
8. Primary display
9. Measuring unit for primary display

ACCESSORIES

HI 1296D	Combination preamplified pH/temperature probe with titanium body, DIN connector and 1 m (3.3') cable
HI 7004L	pH 4.01 buffer solution, 500 mL
HI 7006L	pH 6.86 buffer solution, 500 mL
HI 7007L	pH 7.01 buffer solution, 500 mL
HI 7009L	pH 9.18 buffer solution, 500 mL
HI 7010L	pH 10.01 buffer solution, 500 mL
HI 70300L	pH electrode Storage solution, 500 mL
HI 7061L	pH electrode Cleaning solution, 500 mL
HI 77400P	Calibration kit (pH 4 & 7, 20 mL, 5 pcs each)
HI 76405	Electrode holder

SPECIFICATIONS

Range (*)	-2.00 to 16.00 -5.0 to 105.0°C / 23.0 to 221.0°F
Resolution	0.01 pH 0.1°C / 0.1°F
Accuracy (@20°C/68°F)	±0.02 pH ±0.5°C up to 60°C; ±1.0°F up to 140°F; ±0.2°C above 60°C
Temperature Compensation	Automatic -5.0 to 105.0°C (23 to 221°F)
pH Calibration	Automatic choose between 2 points (standard: 4.01/7.01/10.01 or NIST: 4.01/7.01/10.01)
Probe (included)	amplified pH probe with titanium body, DIN connector and 1 m (3.3') cable
Battery Type	AAA
Battery Life	Approximately 1200 hours of use
Auto-off	After 8 min of no use
Environment	0 to 50°C (32 to 122°F); 100% RH
Dimensions	152x 58 x 30mm (6.0 x 2.3 x 1.2")
Weight	210g (7.4 oz.)

(*) the temperature range is limited to 80°C (176°F) the HI 1296D sensor.

To clean the meter, use water only.

Hanna Instruments reserves the right to modify construction and appearance of its products without notice.

IST991

sales@GlobalTestSupply.com

www.GlobalTestSupply.com

Find Quality Products Online at:

ign, nce

1/11

OPERATIONAL GUIDE

First time you use the instrument open the battery compartment and put the batteries inside observing the polarity.

To connect the probe

With the meter turned off, connect the **HI 1296D** probe to the DIN socket on the bottom of the meter by aligning the pins and pushing in the plug. Tighten the nut to ensure a good connection. Remove the protective cap from the probe before taking any measurements.

To turn the meter ON and check the battery status

Press the **ON/OFF/MODE** button until the display lights up. At start-up, all the LCD segments are displayed for 1 second, then the percent indication of the remaining battery life is displayed for another second (E.g. "100 % BATT"). The meter then enters the normal measuring mode.

Note: If the display needs to be checked, keep the **ON** button pressed while turning the meter on. The meter will display all segments as long as the button is pressed.

To freeze the display

While in measurement mode, press the **SET/HOLD** button, "HOLD" appears on the secondary display and the reading will be frozen on the LCD (E.g. "pH 5.73 HOLD"). Press any button to return to normal mode.

To turn the meter OFF

While in measurement mode, press the **ON/OFF/MODE** button. "OFF" will appear on the secondary display. Release the button.

The meter is provided with an acoustic signal active when a key is pressed.

Note: When the meter detects the absence of probe at its input, the message "Probe not connected" appears on LCD and "---" blinking on LCD lines. When a probe is connected, the "Probe not connected" tag is turned off, and the readings are displayed on LCD lines.

pH MEASUREMENT & CALIBRATION

- Make sure the meter has been calibrated before use.
- If the probe is dry, soak it in **HI 70300** storage solution for 30 minutes to reactivate it.
- Submerge the probe in the sample to be tested while stirring it gently. Wait until the "Not Stable" tag on the LCD is turned off.
- The pH value automatically compensated for temperature is displayed on the primary LCD, while the secondary LCD displays the sample temperature.
- If measurements are taken in different samples successively, rinse the probe tip thoroughly to eliminate cross-contamination. After cleaning, rinse the probe tip with some deionized water and some of the sample to be measured.

pH calibration

For better accuracy, frequent calibration of the pH sensor with the meter is recommended. In addition, the meter must be recalibrated whenever:

- a) The pH electrode is replaced.
 - b) After testing aggressive chemicals.
 - c) Where high accuracy is required.
 - d) At least once a month.
- While in pH measurement mode, press and hold the **ON/OFF/MODE** button until "OFF" on the secondary display is replaced by "CAL". Release the button. Place the sensor into the first calibration buffer.
 - The meter enters the calibration mode, displaying "pH 7.01 USE" (or "pH 6.86 USE" if the NIST buffer set was selected).
 - For a single-point calibration place the probe in any buffer from the selected buffer set (E.g. pH 4.01, 9.18 or 10.01). The meter will automatically recognize the buffer value.
 - If the buffer is not recognized or the calibration offset is out of the accepted range "---WRONG" is displayed.
 - If the buffer is recognised "REC" is displayed until the reading is stable and the calibration is accepted.
 - If using pH 7.01 (or pH 6.86 from NIST buffer set), after acceptance of the buffer press any key to exit. "OK1" message is displayed and meter returns to pH measurement mode.
 - If using 4.01 or 10.01 (9.18) buffer the "OK1" message is displayed and meter returns to pH measurement mode.

- For a two-point calibration, proceed as for the single point calibration selecting 7.01/6.86 pH buffer.
- The "pH 4.01 USE" message is displayed.
- Place the probe in the second calibration buffer (pH 4.01 or 10.01, or, if using NIST, pH 4.01 or 9.18). When the second buffer is accepted, the LCD will display "OK2" for 1 second and the meter will return to the normal measurement mode.
- If the buffer is not recognised or the slope is out of accepted range "---WRONG" is displayed. Change the buffer, clean the electrode or press any key to exit calibration.

Note: When the calibration procedure is completed, the "Calibrated" tag is turned on.

It is always recommended to carry out a two-point calibration for better accuracy.

To exit calibration and reset default values

- After entering the calibration mode and before the first point is accepted, it is possible to quit the procedure and return to the last calibration data by pressing the **ON/OFF/MODE** button. The LCD displays "---ESC" for 1 second and the meter returns to normal mode.
- To reset the default values and clear a previous calibration, press the **SET/HOLD** button after entering the calibration mode and before the first point is accepted. The LCD displays "---CLR" for 1 second, the meter resets to the default calibration and the "Calibrated" tag on the LCD disappears.

METER SETUP

Setup mode allows the selection of the temperature unit and the type of pH buffer set.

To enter the Setup mode, press and hold the **MODE** button until "CAL" on the secondary display is replaced by "TEMP" and the current temperature "TEMP °C"). Then:

- for °C/°F selection, use the **SET/HOLD** button. If the temperature unit has been selected, press the **MODE** button to enter the buffer set selection. Press **ON/OFF/MODE** twice to return to the normal mode.
- to change the type of calibration buffer set, after temperature unit, the meter will show the buffer set: "pH 7.01 BUFF" (for standard buffer set 10.01) or "pH 6.86 BUFF" (for NIST buffer set 9.18). Change the set with the **SET/HOLD** button. Press **ON/OFF/MODE** to return to normal mode.

BATTERY REPLACEMENT

The meter is supplied with batteries.

The meter displays the remaining battery percentage on the secondary display. When the level is below 5%, the "BATT" symbol blinks to indicate a low battery condition. If the battery is not replaced soon enough to cause erroneous readings, the Battery End of Life System (BEPLS) turns the meter off. It is recommended to replace the batteries as soon as the display will flash the "BATT" symbol.

To replace the batteries, follow the next steps:

- Open the battery compartment cap (on the back of the instrument).
- Remove old batteries.
- Replace new batteries observing the polarity indicated on the instrument.
- Close the battery compartment cap.